



There is a lamp on the control box to indicate whether the cushion is 'on' or 'off' (green or red) and a push button to turn the cushion on/off.

Check whether the system is on or off

- Press the push button briefly
 - The lamp indicates whether the system is on or off. Red is off, green is on.

Switching on and off

- Hold the push button for 3 seconds.
 - The lamp changes colour and the system is now on or off.
 - This is also indicated by an acoustic signal: 2 short beeps for on, 1 beep for off.
 - When someone sits down on the cushion, the lamp on the control box will show the status of the system.

Reminder function

- When the system is switched off and a client is helped out of the chair, the next time the chair is used the cushion will emit 9 short beeps to remind the caregiver that the system is switched off.
 - This reminder will only be given once, so that the chair can be used, for example by visitors, without activating the system.

By default, the Optiseat works as follows

- Switch the system on.
 - When the client gets out of the chair, a beep sounds for 4 seconds (adjustable), to activate the client's self-correcting ability.
 - If the client does not sit again within 4 seconds, a call signal will be sent to the receiver or the pager.
 - If the client does sit down within 4 seconds, the beep signal stops and nothing else happens.

Batteries

In the Optiseat only two AA 1.5V (penlight) batteries may be used.

Discharged battery notification

When the batteries are discharged, the lamp will flash red/green every 5 seconds. If the system is switched on, a beep will also be produced once every minute.

Changes to the default settings

The technical customer service representative can change the default time settings and the length of the delay before transmission. See the technical manual.

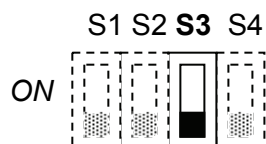
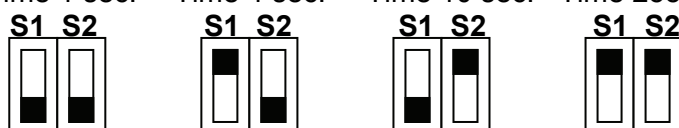
**Opening the unit**

To open the unit you must remove four screws. The two lower screws are located in the battery compartment.

Internal settings**S1 + S2**

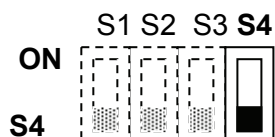
The Optiseat has an alarm delay. This is the time that the person has to sit again before the alarm will sound. Four time settings are available. These times can be changed upon request. The standard delay times you can choose from are: 1, 4, 10 and 255 seconds (other lengths of time are optional). The selection from among these 4 lengths of time is made using the first 2 switches **S1** and **S2**.

Time 1 sec. Time 4 sec. Time 10 sec. Time 255 sec.:

**S3**

The third switch (**S3**) is used to switch the buzzer on and off. If this is up (**ON**), the buzzer is audible until the end of the set alarm delay time.

If **S3** is down (**OFF**), the beep will only sound for one second when the client gets off the cushion, regardless of the alarm delay time.

**S4****S4 is not applicable here.****Discharged battery notification**

When the led flashes red/green every 5 seconds, the batteries are due for replacement. If the system is on, a short beep will also be audible once every minute. Replace the two penlight batteries. If you wait too long before replacing the batteries, a discharged battery notification will appear on the receiver or pager. When multiple Optiseat® cushions are used with a single receiver, the receiver does not indicate which cushion is producing the signal. In this case you must check all the Optiseat® cushions.

When closing the unit, make sure the lamp is positioned in the corresponding opening.